



PRIORITY DOCUMENT

SUBMITTED OR TRANSMITTED IN COMPLIANCE WITH RULE 17.1(a) OR (b) PCI/GB 2004 / D 0 2 8 0 2



INVESTOR IN PEOPLE

The Patent Office Concept House Cardiff Road Newport

South Wales

NP10 8

REG 15 JUL 2004

WIPO PCT

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

I also certify that the attached copy of the request for grant of a Patent (Form 1/77) bears an amendment, effected by this office, following a request by the applicant and agreed to by the Comptroller-General.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed Mag

Dated 8 July 2004

P (atents Form 1/77 THE PATENT OFFICE A 16) -2 JUL 2003	Patent 02JUL03 819 Office P01/770000.0	517-1 C22964 C4031546839
(C)	ee the notes on the back of this form. You can also get an planatory leaflet from the Patent Office to help you fill in is form) Your reference	⊆2 JUL 2003	The Patent Office Cardiff Road Newport South Wales NP10 8QQ
		IBST Ø91	•
_	Patent application number (The Patent Office will fill in this part)	0315468.9	
4.	BEDFORDSHIRE, MK43 CAI I BOULEVARD FLEMING, 2 Patents ADP number (if you know it) If the applicant is a corporate body, give the country/state of its incorporation Title of the invention IMMUNOC FOR MEASUREMENT OF F	CRANFIELD UNIVERSITY, CL, UK AND DIACLONE, P 5020 CEDEX, BESANCON 06458848007 08665986001 MEMICAL METHODS AND FUEL OXYGENATES	o Box 1985, France
5.	Name of your agent (if you have one) "Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)	EMBA TICHELL CRANTICLD UNIVERSITY SUSOE BUDFORDSHIRE MKNS HDT	MENBACH ELLS LES YORK HOUSE, 23 KINGSWAY, LOMPON, PFSI WCZB GMS. 177
 3.	Patents ADP number (If you know it) If you are declaring priority from	07656473002 .	109006 -
	If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (if you know it) the or each application number	Country Priority application number (If you know it)	Date of filing (day / month / year)
	If this application is divided or otherwise derived from an earlier UK application, give the number and the filling date of the earlier application	Number of earlier application	Date of filing (day / month / year)

75S

8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (Answer Yes' 16:

c) any named applicant is a corporate body.

applicant, or

See note (d))

a) any applicant named in part 3 is not an inventor, orb) there is an inventor who is not named as an

Patents Form 1/77

Enter the number of sheets for any of the following items you are filing with this form. Do not count copies of the same document	
Continuation sheets of this form Description	15. CP
Claim(s)	
Abstract	
Drawing(s)	2 only
10. If you are also filing any of the following, state how many against each item.	
Priority documents	
Translations of priority documents	
Statement of inventorship and right to grant of a patent (Patents Form 7/77)	
Request for preliminary examination and search (Patents Form 9/77)	
Request for substantive examination (Patents Form 10/7)	
Any other document (please specify	
11.	I/We request the grant of a patent on the basis of this application.
	Signature L. 3. Give Date 30-06-03
12. Name and daytime telephone number of	LINDA TICWELL
person to contact in the United Kingdom	01525 863168

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

- a) If you need help to fill in this form or you have any questions, please contact the Patent Office on 08459 500505.
- b) Write your answers in capital letters using black ink or you may type them.
- c) If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- d) If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- e) Once you have filled in the form you must remember to sign and date it.
- f) For details of the fee and ways to pay please contact the Patent Office.

Field of the Invention

The present invention relates to immunochemical assays and associated immunoreagent production methods for determination of members of the *tert*-butyl group of fuel oxygenates and associated compounds. The approach, comprising an immunochemical-based assay and associated reagents may be used in either the laboratory of the field for the rapid and routine diagnosis of samples suspected of containing said compounds thereby allowing rapid analytical and diagnostic information to be obtained.

Background to the invention

Human endeavour is hugely dependent upon the consumption of fuels for power generation purposes. For example, road transport is an integral feature of today's society, essential for industry, commerce and recreational activities. Global demand for petroleum necessitates usage of large fuel storage tanks and underground oil pipelines, which present the potential for significant environmental contamination from persistent and hazardous petroleum hydrocarbons. Particularly with underground systems, leaks can remain undetected and can be transported through soils and water systems via natural dispersion mechanisms.

From the early 1970's, increasing concern about public health and environmental impacts from exhaust emissions has led to changes in the formulations of motor fuels. Revised European and US regulations have been targeted to reduce emissions from motor vehicles and improve air quality. Oxygenated compounds have been incorporated into fuels to aid combustion and reduce exhaust emissions responsible for air pollution. These provided a replacement for benzene and have aided fuel

producers in reaching octane levels required by new fuel standards.¹ The choice of oxygenate varies between countries and US states, influenced by the availability of raw materials, individual economies and political status.

Currently, the most commonly used fuel oxygenates are methyl *tert*-butyl ether (MTBE) and ethanol. Other oxygenates in use include ethyl *tert*-butyl ether (ETBE), methyl *tert*-amyl ether (TAME), isopropyl ether, and *tert*-butyl alcohol (TBA). These compounds have been primarily used to achieve mandatory fuel standards throughout the US and Europe. Whilst petroleum may-comprise of >120 different compounds, leaded, unleaded and premium grade fuels may contain >15% oxygenates by volume.

MTBE has since become one of the most significant environmental pollutants in recent years. Sampling results from the US Geological Survey's (USGS) National Air and Water Quality Assessment (NAWQA) Programme during 1990-1998 found high incidences of MTBE contamination in both confined and unconfined ground water aquifers and in drinking water.² Approximately 60% of US drinking water is extracted from surface water systems. Whilst the potential toxicity is still debated, MTBE has an extremely low taste threshold and low-level contamination has led to losses of drinking water supplies. 71% of the water supply of Santa Monica, California, has been tainted with MTBE, requiring the importation of water supplies from outside of the city limits at an annual cost of \$3.5 million per annum.³ MTBE has been

¹ Environment Agency (2001) National groundwater and Contaminated Land Centre Annual Report 1999/2000

² Moran, M.J. et al. (1999) Relations between the detection of MTBE in surface and groundwater and its content in gasoline. US Geological Survey Report, 1999.

³ Environment Agency. The fuel additive MTBE – a groundwater protection issue? National Groundwater and Contaminated Land Centre.

The main source of environmental oxygenate contamination is from oxygenate-blended petroleum, through leaking underground storage tanks (LUSTs), transfer spillage, petroleum facilities and/or accidental spills. Poor management of stored fuel and subsequent leakages have resulted in major, global environmental contamination from petroleum and related components, typically BTEX compounds (benzene, toluene, ethylbenzene and *m-*, *o-*, *p-*xylenes), low weight alkanes, such as *n-*pentane, *n-*hexane, and fuel oxygenates. MTBE has particular environmental significance due to its high solubility in water, with numerous incidences of ground water and/or drinking water contamination. Surface water systems may become contaminated through ground water, atmospheric deposition, storm water run-off and direct releases by industrial and recreational activities.

State of the art

Traditionally, environmental analysis of petroleum and related compounds are performed using laboratory-based methods such as gas or liquid chromatography allied to suitable detection methodologies. The two most widely used analytical methods for detection of fuel oxygenates are EPA Method #8260 (Volatile organic compounds by Gas Chromatography/Mass Spectroscopy (GS/MS) and Method #8015 (Non-halogenated organics using Gas Chromatography/Flame Ionisation Detector (GC/FID). Additionally, there are problems with sample collection methods and sample preservation for these laboratory-based techniques. Whilst a wide range of

⁴ Environment Agency (2001). Fuelling Stations Construction and Operation. Environmental Agency Pollution Prevention Guidelines PPG7.

⁵ Bradley, P.M. et al. (1999). Aerobic mineralization of MTBE and tert-butyl alcohol by stream bed microorganisms. Environmental Science and Technology, 33, pp. 1877-1879.

analytical techniques are available, these remain primarily laboratory-based. Due to escalating reports of petroleum and specifically MTBE contamination across the US and Europe, there remains an urgent need for simple and accurate diagnostic tools for *in-situ* analysis and continuous monitoring applications.

Immunoassays

The specific, sensitive and low-cost decentralised determination of a number of different types of common environmental contaminant has been achieved through the use of immunoassays. Immunoassays allow the detection and measurement of target compounds using specific binding characteristics of antibodies and can be found in a wide variety of formats and are increasingly being developed and employed for environmental monitoring purposes. These have acquired wide acceptance in the USA, with the US EPA recognising and releasing official ELISA methods for the determination of certain compounds, e.g. certain pesticides, PAHs and PCBs. 'Rapid' immunoassay test kits are commercially available for a wide range of analytes, which provide a relatively inexpensive, rapid (<2 h), sensitive screening method for analyte detection, commonly in aqueous and soil matrices. Castillo et al.6 have evaluated the sensitive determination of for the accurate and immunoassay test kits pentachlorophenol, carcinogenic PAHs and BTEX compounds within industrial effluents.

Fundamental principles of immunoassays

The fundamental principle of immunoassays (IAs) is that they utilise biologically generated immunoglobulin proteins – antibodies (Ab) - which react with specific target compounds – namely the target analyte, commonly referred to as the 'antigen'

⁶ Castillo, M. et al. (1998). Environmental Science and Technology, Vol. 32(14), pp. 2180-2184.

(Ag) to form antibody-antigen complexes (Ab-Ag). All IAs are based on the selectivity and sensitivity of this Ab-Ag reaction. Due to the wide range of variations in assay design, there is no universal, consistent classification system, or terminology. Each immunoassay will need to be specifically developed to the target of interest, i.e. the antigen. The Ab is the key reagent and these generally determine the assay's characteristics. Ab-Ag binding arises from structural complimentarity between the 2 molecules, stabilised by binding through a combination of Van der Waals forces, electrostatic interactions, hydrogen bonding and hydrophobic interactions.

Antibody production

The initial stimulation of Immunoglobulin G (IgG) antibodies is achieved through injecting animals, commonly mice, rabbits or sheep, with the target immunogen (antigen). The *in vivo* administration of an immunogen stimulates B-lymphocyte cells to produce and secrete antibodies into the blood stream that are capable of binding to, and with the help of other factors in the immune response, destroying the invading blood-borne foreign body. Since each stimulated B-lymphocyte cell will produce a unique antibody 'clone' exhibiting a specific binding reaction, and hence affinity, to the target analyte, a 'polyclonal' mixture of antibodies derived from all of the individual Ab secreting cells, is elaborated within the serum. This polyclonal antiserum can be utilised at this stage, although the IgG fraction will contain many antibodies of differing specificities, many of which may be irrelevant for IA purposes.⁷ A purification step is generally necessary to increase performance and limit the possibility of activity from those irrelevant fractions.

⁷ Edwards, R. (1996) (Ed.) Immunoassays. Essential data Series. John Wiley and Sons, Chichester, UK.

batch variations in the quality of the antiserum and the source of the preparation will cease on death of the animal host. This is a particular problem in immunoassay manufacture when product consistency is paramount to maintaining a valid analytical/diagnostic tool. This issue has been overcome through the development of monoclonal antibody production methods, which exploit the use of neoplastic multiple myeloma 'tumour' cells — essentially Ab secreting cells that undergo

A problem with polyclonal Ab preparations are that there are significant batch-to-

uncontrolled and rapid cell division. Individual Ab secreting cells are isolated from the animal host, and fused with these myeloma cells to form hybridoma cells that, with careful cultivation, act as effectively immortal cell lines for the production of individual cloned (monoclonal) Ab preparations. Generally, vast numbers of hybridoma cells are generated and then screened for Ab-Ag binding efficacy. The high purity, homogeneity and cloning ability of MAbs enables easier purification and subsequent labelling of these highly specific antibodies. Although costs and practical investments are initially high, these are now the preferred and established practice for IA design.

A significant issue arises in that most organic pollutants are of insufficient molecular weight to engender an immunogenic response. Molecules with a molecular weight of <3000 are not immunogenic and those <5000 may be too weak to induce an adequate immunogenic response. In these cases, the molecules must be conjugated to a much larger carrier protein in order to provoke an immune response. Both MTBE and ETBE have molecular weights of 88 and 102 respectively and are termed 'haptens' (i.e. MW around 100 or less) and thus need to be conjugated in order to stimulate host antibody production. Only those antibodies binding specifically to the haptenic determinant, as opposed to the carrier protein will be of diagnostic use. Theoretically,

various immunogenic carriers can be used, the most common being bovine serum albumin (BSA), human serum albumin (HSA), rabbit thyroglobin and keyhole limpet haemacyanin (KLH).

Immunoassay Design

There are many possible immunoassay formats reported in the literature that exploit the fundamental principle of quantifying the extent of Ab-Ag binding. The method most used for field-based determinative purposes is termed the indirect competitive assay format. The method is indirect in that one of the immunoreagents must be labelled in order to visualise the extent of Ab-Ag binding and hence quantify the amount of target analyte in the sample solution. The choice of label remains largely dependent on the specific characteristics and expected concentrations of the target analyte.

Suitable labels include enzymes, fluorochromes and radioisotopes, the latter being less widely used in recent years. Enzymes can also be selected to convert non-fluorescent substrates to fluorescent products, e.g. alkaline phosphatase. Photoluminescent compounds in immunoassays (FIAs) can provide even greater sensitivity that colorimetric substrates. Fluoroscein isothiocyanate (FITC) is usually the label of choice for immunofluorescent IAs, as coupling procedures are straightforward works with almost all antibodies. Additional labels commonly used include tetramethylrhodamine isothiocyanate (TRITC), Texas Red (TR) or phycoerythrin (PE).

Enzyme Linked Immunoassays (ELISA)

Enzyme linked immunosorbent assays (ELISAs) are based on the combination of selective antibodies with sensitive enzymes that react with a substrate to produce a

detectable colour change, e.g. commonly horseradish peroxidase (HRP), or alkaline phosphatase. Specific enzymes are incorporated to link to target contaminants and through enzymatic actions on the colouring agent (chromogen), enabling both qualitative and quantitative analysis through the catalytic capability of the enzymes. Many regulated contaminants, such as pesticides, polyaromatic hydrocarbons (PAHs) and other organic pollutants can now be detected on-site using available test kits utilising enzyme linked immunoassay technology. A number of solid supports can be used as a means of separation in an ELISA system. Solid support within ELISA systems can be the traditional microtitre plate systems, coated tubes or using covalently bound antibodies to magnetic or latex particles. Due to their relative simplicity, flexibility, speed and cost, IAs are now generally considered an effective and suitable form of analysis, suitable for both laboratory and field diagnostics. Such IA systems are capable of detecting very low levels of contamination, in some cases as low as the ng/L or ppb level.⁸

For petroleum contaminants, a number of IA based field kits are commercially available for on-site analysis, e.g. the Strategic Diagnostics Incorporated RaPID Assay Petroleum Fuel Kits. These field kits utilise Abs bound onto microscopic particles, for which results are interpreted using a portable microprocessor controlled spectrophotometer or photometer. Detection levels can be at µg/L (ppm) or ng/L (ppb) concentrations, depending on the analyte in question and product used. Reportedly, results can be obtained in less than two hours. However, whilst the IA aspects can be accomplished in some cases <60 minutes. A further extraction procedures are necessary for soil analysis, which increase the complexity and time scale of analysis.

Summary of the invention

There is extreme current interest in MTBE and related fuel oxygenates -ETBE, TAME, TBA etc. - hereafter referred to as the *tert*-butyl oxygenate family, and their environmental impact. Despite this fact there is no prior art concerning the development of immunochemical methods for determination of such compounds. It is worthy of note that commercially available test kits to the other major classes of petroleum species, such as the BTEX compounds, TPHs (total petroleum hydrocarbons) and PAHs (polynuclear aromatic hydrocarbons) do currently exist.

This invention therefore pertains to the processes, technologies and associated knowledge associated with the generation of immunoreagents and the development of immunoassays for the antibody-based determination of members of the *tert*-butyl oxygenate family. The benefits of such an approach lie in the simplicity of the assay procedure and amenability of said procedure to decentralised usage. The method offers the following benefits over current MTBE determination methods:

- The simplicity of the procedure and requirement for simple instrumentation (optical reader, electrochemical monitoring device or such other signal interrogation device) renders the method amenable to decentralised operation.
- High sample throughput: achieved through the ability to assay a large number of samples simultaneously.

⁸ http://www.sdix.com

- A high degree of assay specificity and sensitivity due to the binding complimentarity and binding affinity between the antibody and target analyte(s).
- Assay rapidity, typically in the region of 1-2h.
- In situ sample determination removes the requirement for transportation of samples to a centralised facility and the concomitant degradation of sample during transportation and storage.
- Despite the costs associated with antibody preparation, particularly
 monoclonal antibody preparation, the small amounts of Ab material and other
 reagents required for each assay measurement results in low assay costs
 relative to the costs associated with running complex instrumentation at
 dedicated centralised analytical facilities.
- Operator costs are reduced due to the speed/throughput of the assay and the lower levels of training required relative to their laboratory counterpart.

In order to generate antibodies with binding specificities directed towards the *tert*-butyl oxygenate family, account must be made of the extreme low molecular weight of these compounds. Research in our laboratories has focussed on producing *tert*-butyl oxygenate-protein conjugates capable of engendering an immunological response in a host animal species.

oxygenate analogues for immunogen production. These analogues may be polymers (dimers, trimers, polymers etc.), consisting of polymerisation of two or more tertbutyl oxygenate monomers, to form extended repeated tert-butyl oxygenate polymer chains, or alternatively, individual molecules exhibiting structural and functional characteristics of the tert-butyl ether target compound in association with a second entity that acts as a spacer between the tert-butyl ether functionality and associated carrier compound. Additionally, these compounds will also contain a functionality located away from the tert-butyl ether moiety that can be used for carrier compound conjugation purposes, such as conjugation to suitable carrier proteins. Such entities, henceforth referred to as [tert-butyl oxygenate]n-carrier protein conjugates (i.e. tertbutyl ether polymer structures and tert-butyl ether-spacer compounds) were found to elicit the desired immunological response in host animals, with recovery of antibody with specificity and high binding affinity for both the tert-butyl oxygenate polymer and original monomer compound. Simple conjugation of members of the tert-butyl oxygenate family to standard carrier proteins yielded antibodies with significantly weaker antigen binding affinity.

This invention further relates to the use of methodologies whereby the specific binding capabilities of these antibodies, produced using either the *tert*-butyl oxygenate polymer-carrier protein or *tert*-butyl oxygenate monomer-carrier protein immunogen route, are exploited for the creation of immunodiagnostic methods for the determination of members of the *tert*-butyl oxygenate family. The immunoassay format may vary in nature but relates to the use of immunochemical assay formats and specific assay labels that are able to visualise the extent of Ab-Ag binding and hence lead to the determination of the *tert*-butyl oxygenate family members.

The assay format may be competitive or non-competitive in operation and may include alternative embodiments of either approach, such as the use of particulates, such as latex beads, magnetic beads and the well described lateral flow assay format. The assay label may be linked directly to one or more of the immunoassay reagents, such as antibody or antigen analogue, or may be introduced by alternative means, such as via the well-known streptavidin-biotin binding complex, or through conjugation of the label to a second antibody or binding component with binding specificity directed towards structures on the primary anti-tert-butyl oxygenate antibody. The label may be, but is not limited to an enzyme, chromophore, fluorophore or other optically active agent, chemically active agent, electrochemically active agent or other such suitable compound in which the specific properties of the label can be used to visualise the Ab-Ag binding process. In the case of enzyme labels, the consumption of active substrate and/or the generation of active product or other enzyme mediated effect may be used to produce the assay response. Selection of the transduction process which will be dependent upon the selected assay label may include, but not be limited to optical, electrochemical or other appropriate method. Integration of the immunoassay with appropriate transduction methodologies to produce dedicated sensing tools is a further clear embodiment of the invention described herein.

Whilst it is evident that the invention is suitable for decentralised determination of members of the *tert*-butyl oxygenate family using the indirect immunochemical assay format, it is evident that the methodologies herein described may be equally applied to direct immunochemical assay methods. Direct methods are able to visualise the Ab-Ag binding process without the aid of associated assay labels by use of appropriate transduction methodologies. Such methodologies would include, but would not be

limited to: surface plasmon resonance (SPR) devices, evanescent wave devices, quartz crystalline microbalance (QCM) devices, surface/bulk acoustic wave devices, field-effect transistors or any other such methodology that may be considered direct in operation.

Specific example

A specific embodiment of the invention is now given by way of the following example.

An indirect competitive immunoassay with specificity towards MTBE has been constructed and tested with representative samples. The assay is centred on the competition between microtitre well wall immobilised MTBE-spacer-BSA polymer (7-methoxy-3,7-dimethyloctanal-BSA - the 'coating antigen'), free MTBE in the sample and free anti-MTBE/MTBE-spacer antibody.

Anti-MTBE/MTBE-spacer antibody was first generated by conjugating 7-methoxy-3,7-dimethyloctanal to BSA carrier protein. The resultant conjugate was injected into a mouse host and antibody synthesising cells harvested at the appropriate time. The antibody synthesising cells were fused with myeloma cells to form hybridoma cells which were cultivated and the antibody specificity of each hybridoma assessed with respect to binding to MTBE and 7-methoxy-3,7-dimethyloctanal. Those hybridomas expressing antibody of the required specificity were propagated and used as a source of anti-MTBE/7-methoxy-3,7-dimethyloctanal antibody for immunoassay applications.

A competitive immunoassay is described in this particular embodiment. Free MTBE antigen present in the sample was allowed to compete with the immobilised coating

antigen for anti-MTBE/7-methoxy-3,7-dimethyloctanal antibody. The quantitative

the assay is dependent upon free MTBE-Ab binding. In essence, the lower the concentration of free MTBE in the sample solution, the greater the extent of antibody binding to the coating antigen. On washing the microtitre wells, free MTBE-antibody complexes are displaced from the well. An antibody-specific enzyme tracer is then added to the system. Following the binding and washing steps, residual enzyme activity, which is inversely proportional to the original free MTBE concentration, is then determined by adding enzyme substrate and recording the optical density of the coloured enzymic product after a suitable period of colour development. Low concentrations of MTBE in the sample gave an inversely proportional large signal output, whilst high concentrations of free MTBE in solution resulted in a reduced assay signal.

Microtitre well preparation

Wells were prepared by adding 100µl of coating antigen (500ng/ml) and then incubated at 4°C for 24 hours. Plates were then aspirated and washed with wash solution (0.01% v/v Tween 20 in Reverse Osmosis (RO) water) at 300µl/well. After each wash step, plates were inverted and blotted against clean paper towelling. Unless otherwise stated, plates were blocked using 5% v/v BSA blocker buffer and incubated at room temperature (RT) for 2 hours. Plates were then aspirated and two wash steps performed. Dried plates were then sealed and stored with desiccant at 4-8°C prior to use.

ELISA immunoassay procedure

A 100µl volume of MTBE-containing sample was added to each microplate well followed by 100µl of biotinylated anti-MTBE/7-methoxy-3,7-dimethyloctanal Ab.

step, as described for plate preparation, was repeated $\times 4$. Streptavidin-HRP (100 μ l per well) was added, the plates sealed and incubated at RT for a further 30 minutes. Aspiration and 4 further wash steps were made. To each well, 100 μ l of the HRP enzyme substrate TMB (3,3',5,5' tetramethyl benzidine) was added. The plates were then covered and incubated for 10 minutes. At the end of the colour development time, the colour reaction was stopped by the addition of 100 μ l 2M H₂SO₄ and the OD₄₅₀ immediately measured. Any necessary dilutions were made with RO water.

Results

Optimal concentrations of both Ag and Ab (biotinylated and native) were determined by checkerboard assay (Figure 1). The optimal concentrations of coating antigen conjugate and biotinylated Ab were all found to be 500ng/ml.

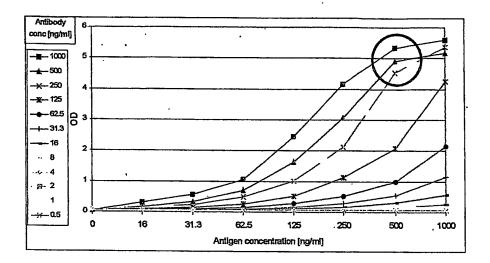


Figure 1. Checkerboard assay for the determination of optimal coating Ag and primary Ab concentrations. The highest signal was achieved using 500ng/ml coating antigen and 500µg/ml anti-MTBE/7-methoxy-3,7-dimethyloctanal antibody.

A 500 ppm MTBE stock solution was prepared in 100 mM phosphate buffered saline (PBS, pH 7.4) and serially diluted (1/10) in 100 mM PBS to a final concentration of 0.5 ppb. A competitive ELISA was then performed. Figure 2 shows a reproducible sigmoidal curve of assay response (OD_{450}) as a function of increasing free MTBE concentration. As the ELISA is a competitive assay, the OD_{450} is inversely proportional to increasing concentrations of free MTBE. The assay was repeated on different days and with fresh solutions of all reagents. Results are shown in Figure 3.

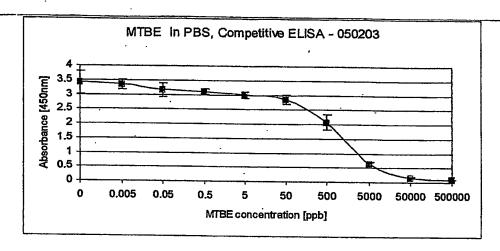


Figure 2. Standard ELISA result comparison with spiked MTBE samples in PBS. The dynamic range of the assay was found to be 50-5000 ppb for free MTBE.

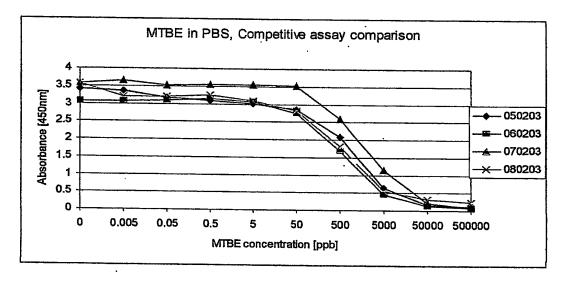


Figure 3. Standard ELISA result comparison with spiked MTBE samples in PBS. The dynamic range of the assay was found to be 50-5000 ppb for free MTBE.. Plates 050203 and 060203 were blocked with 5 % v/v BSA blocking buffer. Plates 070203 and 080203 were blocked with 1% BSA blocking buffer.

This Page is Inserted by IFW Indexing and Scanning Operations and is not part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

□ BLACK BORDERS
□ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
□ FADED TEXT OR DRAWING
□ BLURRED OR ILLEGIBLE TEXT OR DRAWING
□ SKEWED/SLANTED IMAGES
□ COLOR OR BLACK AND WHITE PHOTOGRAPHS
□ GRAY SCALE DOCUMENTS
□ LINES OR MARKS ON ORIGINAL DOCUMENT
□ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY

IMAGES ARE BEST AVAILABLE COPY.

OTHER:

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.